

neoligaments™



## Poly-Tape

### High Strength Sterile Polyester Tape

#### Benefits

No requirement to harvest a soft tissue graft to reconstruct the ligament or tendon: reduces morbidity

Open weave structure acts as a scaffold allowing tissue ingrowth that protects the Poly-Tape from abrasion and adds additional strength as it matures

Large range of sizes available to suit varying patient anatomy

Heat sealed ends prevent fraying and facilitate easier passage through bone tunnels and soft tissue

Made from woven high tenacity polyester providing a soft flexible structure with high strength

## Indications


Poly-Tapes are intended for tissue approximation, including use in reconstructing damaged or torn ligaments, tendons and other soft tissues according to the surgeon's own preferred technique, and at his/her discretion. Particular indications include:

- Quadriceps tendon reconstruction, 30mm x 800mm
- Patellar tendon reconstruction, 30mm x 800mm
- Achilles tendon reconstruction, 20mm x 800mm
- Stabilization of the dislocated peroneal tendon, 10mm x 500mm
- Stabilization of the dislocated acromioclavicular joint with associated disruption of the coracoclavicular ligaments, 20mm x 500mm

## Ordering Information

### Standard Poly-Tape (supplied sterile)

Ref	Weave	Width (mm)	Length (mm)	Strength (N)
102-1056	Open	5	500	300
102-1086	Open	5	800	300
102-1057	Open	7	500	400
102-1087	Open	7	800	400
102-1010	Open	10	500	480
102-1080	Open	10	800	480
102-1011	Open	15	500	630
102-1081	Open	15	800	630
102-1012	Open	20	500	950
102-1082	Open	20	800	950
102-1013	Open	30	500	1380
102-1083	Open	30	800	1380
102-1014	Open	40	500	1850
102-1084	Open	40	800	1850
102-1015	Open	50	500	2250
102-1085	Open	50	800	2250
102-1053	Dense	10	500	740

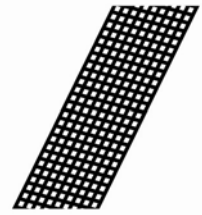

### Tube-Tapes (supplied sterile)

Ref	Weave	Width (mm)	Length (mm)	Strength (N)
102-1040	Open	7 (dia. 5)	500	1000
102-1045	Open	7 (dia. 5)	750	1000
102-1047	Dense	15 (dia. 8)	500	1450



All rights reserved. © Neoligaments 2009. Worldwide patents and patents pending. Neoligaments and Xiros are trademarks of Xiros.



Developed and manufactured by **Neoligaments™**  
A division of Xiros™  
Springfield House Whitehouse Lane Leeds LS19 7UE  
Tel. +44 (0) 113 238 7202 Fax. +44 (0) 113 238 7201  
enquiries@neoligaments.com  
www.neoligaments.com